

Snap-fit Spheres

The Sphere Specifically Designed to Work
with Your Navigation System

The industry's original—and world's leading—passive sphere, the NDI Passive Sphere is a sterile, single-use reflective sphere for use in Image-Guided Surgery (IGS) where it provides **pinpoint accuracy**.

NDI Passive Spheres are packaged to meet the needs of busy OR staff. **Procedure Trays** contain the exact number of spheres required in the most common IGS procedures. Spheres are situated in **deep-seated tray pockets** with a protective lid to minimize the risk of spheres being dropped or popping out when the package is opened. And the spheres are easily picked out of trays – no dumping is required.

The NDI Passive Sphere mounts securely to standard snap-on posts with an audible click for consistent placement. Reliable, accurate, and secure, the NDI Passive Sphere is the reflective passive sphere of choice for **optimal tracking performance**.

Reliable

Consistent sphere shape, surface, and placement are why the NDI Passive Sphere has been the sphere of choice since 1997.

Accurate

The spheres used to accurately and precisely calibrate the IGS System's camera.

Secure

Mounts to standard snap-on posts with an audible click to ensure secure attachment.




Consistent Shape

Consistent Surface

Consistent Placement



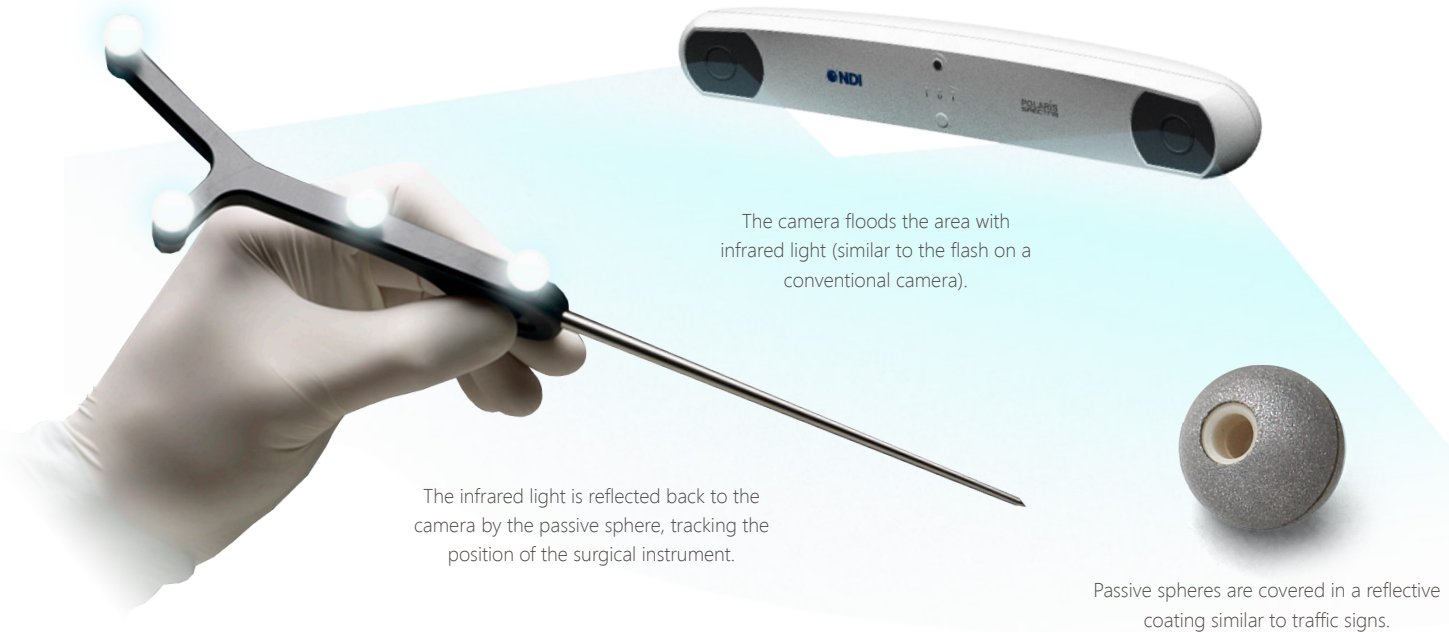

The camera floods the area with infrared light (similar to the flash on a conventional camera).

The infrared light is reflected back to the camera by the passive sphere, tracking the position of the surgical instrument.

Passive spheres are covered in a reflective coating similar to traffic signs.

Procedure Trays



Part Number	Spheres Per Tray	Trays Per Box	Total
8800965	9	2	18
8801320	10	2	20
8800967	12	2	24
8800966	15	2	30

Instrument Trays



Part Number	Spheres Per Tray	Trays Per Box	Total
8801053	3	12	36
8801054	4	12	48
8801002	5	12	60

Pouches



Part Number	Spheres Per Pouch	Pouches Per Box	Total
8800737	1	5	5
8800738	3	5	15
8800739	4	5	20

NDI
spheres.ndigital.com

NDI INTERNATIONAL HEADQUARTERS

CANADA
+ 1 (877) 634-6340
spheres@ndigital.com
www.ndigital.com

NDI EUROPE GMBH

GERMANY
+ 49 (77 32) 8234-0
spheres-eu@ndigital.com
www.ndieurope.com

NDI ASIA PACIFIC

HONG KONG
+ (852) 2802 2205
spheres-ap@ndigital.com
www.ndigital.com